# SUPLITO SUPLITO TWINTOP

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## Table of contents

Preface .....	7
Declaration of conformity .....	7
Service life .....	7
Disposal .....	7
Identification plate .....	7
Reuse.....	7
Warranty .....	7
Parts of the product .....	8
Indication   Contraindication .....	8
Scope of delivery – SuppliTo.....	8
Scope of delivery – SuppliTo twinTop.....	8
Safety instructions.....	8
Assembly – SuppliTo .....	9
Assembly – SuppliTo twinTop .....	9
Adjusting the tilt angle.....	9
Height adjustment of the table top .....	9
Securing the standing position .....	9
Cleaning   Maintenance.....	9
Technical data.....	9

## Preface

Thank you for choosing a product from Drive DeVilbiss.

You will not be disappointed by the design, functionality and quality of this product.

Please read these operating instructions carefully prior to initial use. They contain important safety information and valuable tips for the proper use and care. If you have questions or need additional information, contact your Drive DeVilbiss specialist dealer who supplied your product.

These operating instructions contain all of the important information for adjustment and operation. The most current operating instructions are available as a PDF online or can be requested from customer service.

Enlarged formats for visually impaired persons can be displayed. Blind patients must be instructed by their attendant. Repairs as well as some adjustments require special technical training and must therefore be made by the medical supplies dealer.

## Declaration of conformity

For the product described, DRIVE MEDICAL GMBH & Co. KG declares compliance according to the medical device regulation (EU) 2017/745.

## Service life

Our company assumes a service life of five years for this product as long as the product is used as intended and all maintenance and service requirements are observed.


This service life may be significantly exceeded if the product is carefully handled, repaired, maintained and used and there are no technical limits based on the further development of science and technology.

The service life may also be significantly shortened through extreme and improper use, however. Determination of the service life by our company does not represent any additional guarantee.

## Disposal

If you no longer need the product, please contact your specialist dealer or bring the product to your local disposal center.

## Identification plate

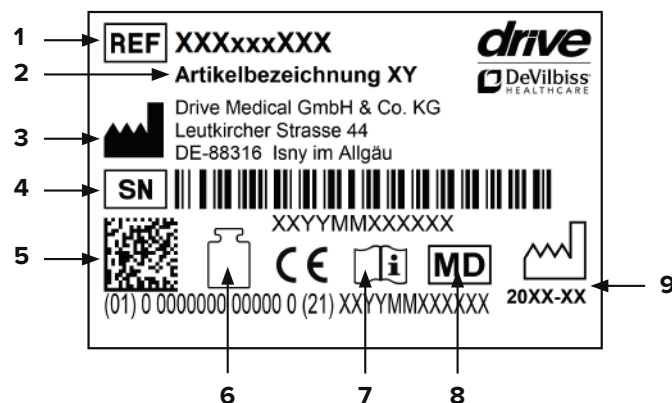1 Article number | 2 Article name | 3 Manufacturer | 4 Serial number | 5 UDI code | 6 Max. load | 7 Follow the instructions for use | 8 Medical device | 9 Production date (year-month)

## Reuse

This product can be reused.

Please note that all technical documents needed for safe handling must also be given to the new user.

The product must be checked by the specialist dealer, cleaned, and restored to optimum condition.

Always follow the manufacturer's instructions during preparation; these are made available on request.

Accessories are available from your specialist dealer or at [www.drivedevilbiss.de](http://www.drivedevilbiss.de)

## Warranty

The warranty period for the product supplied by us is 24 months from the date of purchase. Defective parts will be replaced free of charge if it can be shown that there is a defect in materials or manufacturing.

Parts subject to wear are generally not covered during the normal guarantee period unless the parts need repair or replacement which became necessary as a clear, direct consequence of a manufacturing or material error.

We reserve the right of error and modifications on the basis of technical improvements and changes in design.

### EXCLUSION OF WARRANTY

Parts subject to wear include (among others):  
Casters, brakes

## Parts of the product

### SUPPLITO

1. Table top with connecting tube
2. Table top support
3. 2 x star screw
4. L bar
5. Base with casters

### SUPPLITO TWINTOP

1. TwinTop table top
2. Base frame made of steel tubing (preassembled)
3. L-bar (for assembling the table top 60x40cm)

## Indication | Contraindication

### INDICATION:

The service and side table can be used in combination with a patient or hospital bed. It is intended for users as well as patients and persons in need of care as a valuable aid, for example, during body care, when eating and drinking in bed, or also during therapeutic measures.

### CONTRAINDICATION:

Persons with significantly limited cognitive abilities should use the service and side table only under supervision.

## Scope of delivery – SuppliTo twinTop

1 x	Twintop table top 60 x 40 cm
1 x	Twintop table top 20 x 40 cm
1 x	Base frame made of steel tubing (pre-assembled)
1 x	L bar (for assembling the table top)
3 x	star grip screw
1 x	Instructions for use

## Scope of delivery – SuppliTo

2 x	base with casters
1 x	L bar
1 x	L-bar with table top support
1 x	table top with connecting tube
2 x	star screw
2 x	hexagon socket screw with washer
1 x	spiral spring
1 x	socket head wrench key
1 x	Instructions for use

## Safety instructions

- Before assembly, check the delivery contents, in particular the assembly materials.
- Regularly check the stability and sturdiness of the product and, where applicable, the secure fit of the screws!
- The maximum central load of the service and side table of 10 kg in the horizontal position must be complied with!
- Do not use the service and side table to assist in getting up from a lying position!
- Do not lean or support yourself on the service and side table!
- When stationary, always secure the braking casters.
- Use by children only under adult supervision – risk of pinching and tipping!
- Do not allow hot objects to come into contact with the table top!

## Assembly – SuppliTo

Secure both bases with casters (5) with the L bar (4). Ensure that the feet are mounted at the correctly marked sites (A to A, B to B, C to C and D to D).

Push the table top support (2) into the vertical part of the L bar (4).

Mount the black star screw (3) in the threaded hole provided and tighten the star screw for height adjustment.

Lead the spiral spring over the connecting tube of the table top. Push the table top with connecting tube (1) horizontally onto the table top support (2) and tighten using the star screw (3).

## Assembly – SuppliTo twinTop

The service and side table SuppliTo twinTop is delivered already partially assembled.

Connect the L-bar (3) to the base frame (2), fix it with the star screw. Lead the table top 60 x 40 cm over the L-bar, fix it with the star screw. Push the table top 20 x 40 cm opposite into the L-bar, fix with the star screw.

## Adjusting the tilt angle

### SuppliTo

To change the tilt angle, loosen the star screw with one hand by rotating counterclockwise and with your other hand, adjust the tilt angle as desired. Then tighten the star screw once again by rotating it clockwise.

### SuppliTo twinTop

To change the tilt angle, release the adjustment lever with one hand and adjust the tilt angle with your other hand. The tilt angle can be adjusted in 7 different positions.

## Height adjustment of the table top

### SuppliTo

The height of the table top can be continuously adjusted between 72 and 116 cm. To adjust the height, loosen the **star screw** by rotating it counterclockwise.

Pull the central bar upwards or downwards to the desired height and secure this position by rotating the star screw clockwise.

### SuppliTo twinTop

The height of the table top can be adjusted in 10 positions between 71 and 93 cm. To adjust the height, loosen the star screw by rotating it counterclockwise.

Pull the fixation lever out and adjust the table top to the desired height. Now allow the fixation lever to click into place in the set position. In addition, secure this position by tightening the startscrew clockwise.

## Securing the standing position

Secure the service and side table by using your foot to secure the braking casters by pressing down on the brake.

## Cleaning | Maintenance

The service and side table can be cleaned and disinfected manually. To do so, use a household cleaner and a damp cloth. For disinfection, use a disinfectant listed by the German Society for Hygiene and Microbiology (DGHM)/Association for Applied Hygiene (VAH).

The service and side table does not require any special maintenance.

## Technical data

Product	SuppliTo	SuppliTo twin-Top
white	860000001	860000002
wooden look	860000003	860000004
HMV (medical aids register) no.	50.45.04.0025	50.45.04.0024
Length of table top	60 cm	60 and 40 cm
Width of table top	40 cm	40 and 20 cm
Overall length	62 cm	82 cm
Overall width	40 cm	40 cm
Height adjustment	72 - 116 cm	71 - 93 cm
Tilt angle	0-90° continuous	0-55° in 7 positions
Material	Steel frame, wood table top (plastic-coated)	Steel frame, wood table top (plastic-coated)
Maximum central load on the table top	10 kg	10 kg

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Drive Medical GmbH & Co. KG